

# Assembling Instructions


1- Apply the window self-adhesive seal tape internally to the base frame.



2- Carefully cut the tape in order to perfectly match the opposite end and press well again the tape all around to gain the best adhesion.



3- Take a microfiber vacuum bag (better if a certified Hepa anti-allergen one) and measure two sides far from the vacuum adapter. Draw a rounded line with a pen as in the picture.

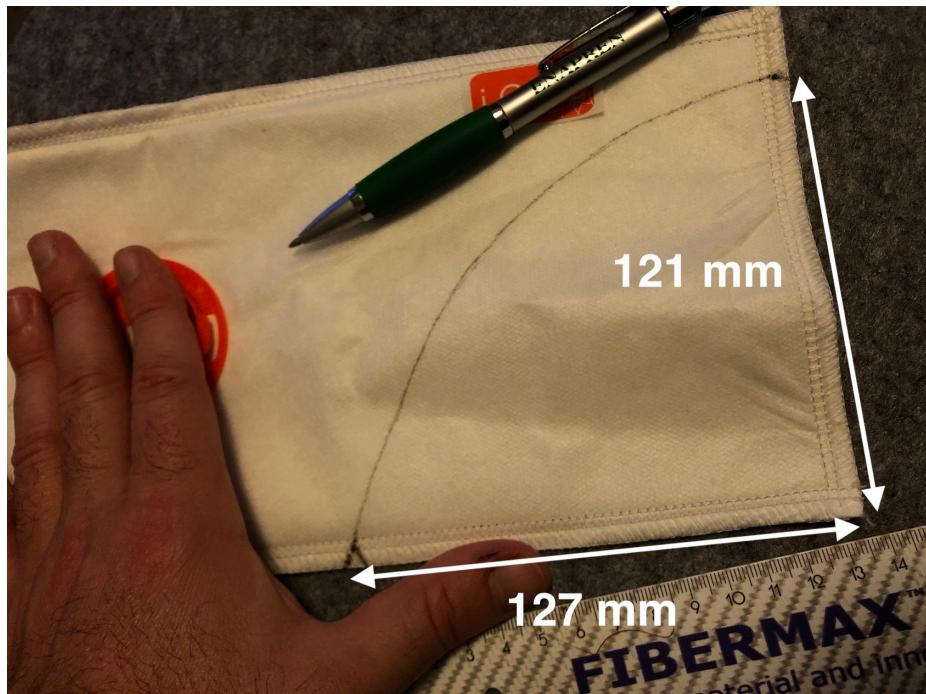

4- Cut the bag following the drawn line.



5- Apply the cone-created filter on the base frame paying attention to perfectly overlap to the frame itself without covering the 3 locking elements. Eventually cut out the material in excess.



6- Lay down from above the top frame onto the filter and start securing the 2 frames together from the lower locking element.



7- Press together the corresponding side-locking elements of base and top frames on both sides.



8- Eventually cut out superfluous filter material. Then squeeze the frames together on the upper side of the mask until the top frame overlaps the base one.



9- Tuck an elastic flat fabric tape (max 1 cm wide) through the first buckle and pass it from inside outward in one slot of the frame and back into the buckle. Repeat this operation for the other end of the elastic with the second buckle.



Remember to sterilize the mask every time after using and I suggest to create more spare filters from vacuum bags.



That's it, I hope it fits on your face.

Mauro.